

Contents

Easycure short-wire ERCP products:

- Adaptor
- Guidewires
- Cannulas
- ProCut Knife
- Sphincterotomes
- Balloons and Baskets
- Stents

Standard long-wire ERCP accessories:

- Cannulas
- Guidewires
- ProCut Knife
- Sphincterotomes
- Balloons and Baskets
- rotaCRUSH Lithotripter
- Endopockets
- Stents
- ENBD Tubes
- Grasping Forceps
- Contrast Syringes





EASYCURE SHORT-WIRE ADAPTOR AND GUIDEWIRE

The new EasyCure ERCP system is for use with short guidewires. The improved wire control ensures a smooth procedure for physicians and nurses alike.

FEATURES

- Continuous fixation of the hydrophilic guidewire during the whole ERCP procedure ensures complete control of the guidewire.
- There is less manipulation at the papilla due to the continuous wire fixation.
- Open system solution – all instruments can be ordered separately and can be used in conjunction with other manufactures short wire systems.

Part number	Description	Pack size
SPE1-Y7	Guide wire locking device for Olympus and Fujinon duodenoscopes	5
WIR1-F9-35-260	gSlider 0.035" hydrophilic guidewire. 260cm long	2

EASYCURE SHORT-WIRE ERCP SYSTEM CANNULA


The new EasyCure ERCP system is for use with short guidewires. The improved wire control ensures a smooth procedure for physicians and nurses alike.

FEATURES

- The guidewire opening channel is split to the distal end potentially providing faster and simpler exchange of instruments.
- Continuous fixation of the guidewire during the whole ERCP procedure ensures complete control of the guidewire.
- The system allows the use of short and long guidewires.



Part number	Description	Pack size
CAN1-B4-18-200-E35	Double lumen cannula tapered end with radiopaque marker	1

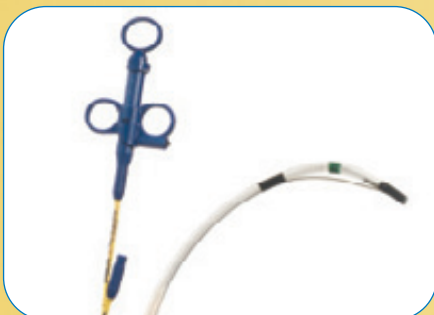

EASYCURE SHORT-WIRE INSULATED PROCUT KNIFE

The new EasyCure ERCP system is for use with short guidewires. The improved wire control ensures a smooth procedure for physicians and nurses alike.

FEATURES

- The guidewire opening channel is split to the distal end.
- Insulated cutting wire provides high current density and therefore reduces pancreatitis and haemorrhage*
- The system allows the use of short and long guidewires.

Part number	Description	Pack size
PAP1-F8-30-E35-OL	ProCut insulated cutting wire EasyCure sphincterotome. Triple lumen, short tapered nose, 30mm monofilament cutting wire (half covered). 0.035" compatible	2

*The Safety Sphincterotome: The device, the technique and preliminary results. DF Martin et al. Endoscopy 1998; 30: 375-378



EASYCURE SHORT-WIRE BALLOON AND BASKET

The new EasyCure ERCP system is for use with short guidewires. The improved wire control ensures a smooth procedure for physicians and nurses alike.

FEATURES

- The guidewire opening channel is split to the distal end.
- Re-sizable balloon and 'Twister basket'.
- Continuous fixation of the guidewire during the whole ERCP procedure ensures complete control of the guidewire.
- The system allows the use of short and long guidewires.

Part number	Description	Pack size
BAS1-H1-20-27-200-E35	Hexagonal 20mm braided wire basket	1
BAL1-C3-18-07-200-E35	Triple lumen biliary extraction balloon catheter. 7fr catheter with 18mm balloon	5

EASYCURE SHORT-WIRE STENT INTRODUCTION KITS

The new EasyCure ERCP system is for use with short guidewires. The improved wire control ensures a smooth procedure for physicians and nurses alike.

FEATURES

- Continuous fixation of the guidewire during the whole ERCP procedure ensures complete control of the guidewire.
- Open system solution - all instruments can be ordered separately and can be used in conjunction with other manufactures short wire systems.



Part number	Description	Pack size
PR01-H7-05-205-E32	One-Step introduction system for 8.5Fr biliary stents	1
PR01-H7-06-205-E35	One-Step introduction system for 10Fr biliary stents	1



EASYCURE SHORT-WIRE SPHINCTEROTOMES

The new EasyCure ERCP system is for use with short guidewires. The improved wire control ensures a smooth procedure for physicians and nurses alike.

FEATURES

- The guidewire opening channel is split to the distal end potentially providing faster and simpler exchange of instruments.
- Continuous fixation of the guidewire during the whole ERCP procedure ensures complete control of the guidewire.

Part number	Description	Pack size
PAP1-C8-30-E35-OL	Triple lumen sphincterotome. Short tapered nose, 30mm monofilament cutting wire	2
PAP1-K8-06-E35-OL	Triple lumen, 6mm needle knife	2



ERCP CANNULAS

The CAN1 Series of ERCP cannulas offer excellent cannulation capabilities covering most clinical situations. All models feature a fluoroscopic tip providing X-ray visibility.

FEATURES

- Single lumen design accepts a 0.035" guide wire for ease of use.
- The pre-curved tip design allows the tip to be curved towards the common bile.
- Visual distal markers aid awareness of depth of cannulation.

Part number	Guide wire	Tip size	Type	Pack size
CAN1-A3-18-200-35	0.035"	5Fr	Standard	10
CAN1-A5-18-200-35	0.035"	5.5Fr	Standard metal	10
CAN1-B3-18-200-35	0.035"	4.3Fr	Tapered	10
CAN1-D5-18-200-35	0.035"	5.5Fr	Bottle tip	10

ERCP gSlider GUIDEWIRES

The gSlider ERCP guidewires are ultra-hydrophilic and have an extremely bright tungsten impregnated 5cm floppy tip. This provides enhanced cannulation and greater fluoroscopic visibility.

FEATURES

- The Nitinol core is kink resistant and thus aids device exchange.
- The gSlider wire incorporates 1cm green and white visual markers to help detect inadvertent wire movement.
- The insulated PTFE jacket enables a sphincterotomy knife to cut over wire without the need to remove it.



Part number	Size (inches)	Length (cm)	Type	Pack size
WIR1-F9-35-260	0.035"	260cm	Standard	2
WIR1-F9-35-450	0.035"	450cm	Standard	2



PROCUT INSULATED SPHINCTEROTOMY KNIFE

The ProCut insulated sphincterotomy knife creates high current density and therefore reduces pancreatitis*, reduces incidence of burn-out on the endoscope and enables cut to be undertaken within a diverticula.

FEATURES

- PTFE coating on the proximal end of the cutting wire provides high current density.
- Precise wire orientation always between 11 and 1 o'clock provides extra safety during sphincterotomy.
- The ultra smooth 4.5fr tapered catheter aids cannulation.

Part number	Nose Length (mm)	Wire Length (mm)	Type	Pack size
PAP1-F8-30-35-OL	6	30	Triple lumen (15mm insulated cutting insulated) wire sphincterotome.	2

*The Safety Sphincterotome: The device, the technique and preliminary results. DF Martin et al. Endoscopy 1998; 30: 375-378



ERCP SPHINCTEROTOMY KNIVES

Cutting takes place at the 11 o'clock position which aids consistent and repeatable cannulation to ensure safer sphincterotomy.

FEATURES

- Precise wire orientation always at the 11 o'clock position provides extra safety during sphincterotomy.
- 4.5fr smooth tapered catheter aids cannulation even for the triple lumen knife.
- Tantalum radiopaque metal marker for clear identification of tip placement.

Part number	Guide Wire Compatibility	Nose Length (mm)	Wire Length (mm)	Type	Pack size
PAP1-C8-20-35-OL	0.035"	6	20	Triple lumen knife	2
PAP1-C8-30-35-OL	0.035"	6	30	Triple lumen knife	2
PAP1-M2-20-35-OL	0.035"	6	20	Double lumen	2

ERCP PRE-CUT SPHINCTEROTOMY KNIVES

Cutting takes place at the 11 o'clock position which aids consistent and repeatable cannulation to ensure safer sphincterotomy.

FEATURES

- Precise catheter orientation always at the 11 o'clock position provides extra safety during sphincterotomy.
- Smooth tapered catheter aids cannulation.
- Tantalum radiopaque metal marker for clear identification of tip placement.



Part number	Guide Wire Compatibility	Nose Length (mm)	Wire Length (mm)	Type	Pack size
PAP1-K1-06-OL	-	-	6	Single lumen needle knife	2
PAP1-K2-06-35-OL	0.035"	-	6	Double lumen needle knife	2



ERCP BALLOON EXTRACTION CATHETER

By combining innovative design features and high quality materials, this product is one of the best balloon extraction catheters available.

FEATURES

- Excellent balloon strength to guarantee efficiency and security of extraction.
- Distal metal markers for fluoroscopic identification of balloon position.
- Coloured graduation markings on distal end give excellent visualisation during video-endoscopy procedures.

Part number	Type	Balloon inflated Ø (mm)	Recommended Guidewire	Pack size
AS1420720	Triple lumen	16	0.035"	5
AS0030520	Double lumen	12	0.025"	5
AS1100720	Double lumen	16	0.035"	5



GRASPING BASKETS

The ERCP BAS1 basket design provides a 'rounder' shape to facilitate opening and stone capture.

FEATURES

- The superior 16 sided enables excellent opening dynamic.
- Baskets incorporate an injection port for contrast injection.
- Built in safety point enabling emergency lithotripsy place if required.
- Unique 'Twister' basket rotates 130° upon closure facilitating stone capture.

Part number	Channel size	Length (cm)	Type and size of Basket	Pack size
BAS1-A2-30-23-200	2.8mm or larger	200	30mm 4 wire hexagonal	2
BAS1-H1-20-27-200-E35	2.8mm or larger	200	30mm Wire guided 4 wire hexagonal	1
BAS1-K2-25-23-200	2.8mm or larger	200	25mm 4 wire spherical twister basket	2

rotaCRUSH ERCP STONE CRUSHING LITHOTRIPTOR


The unique all-in-one design of the rotaCRUSH minimises the complicated assembly process and speeds up the procedure. The thick braided crushing wires enable the basket to withstand a force of up to 300N and the 130° rotation of basket gives superior catching ability.

FEATURES

- All-in-one design incorporating the basket, cannulating sheath, coil sheath and pre-assembled handle.
- Injection port allows for the injection of contrast medium to confirm correct location.
- 'Twist and catch' basket facilitates stone capture.



Part number	Channel size	Type	Pack size
LIT3-K2-30-35-200	4.2mm	All-in-one slide type with inner teflon catheter	1

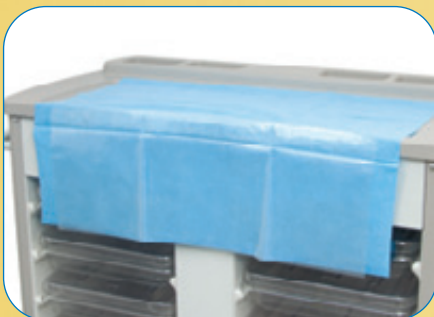

ENDOPOCKETS ERCP TROLLEY ACCESSORY DRAPE

The new EndoPockets not only provide a safe storage place for your different accessories but also gives a sterile cover for your ERCP trolley.

FEATURES

- Available in both 3 and 4 clear plastic pockets for your accessories.
- Triple layer water resistant fabric.
- Fixes to the trolley by a peelable adhesive patch.

Part number	Description	Number of pockets	Pack size
MIS1-Towel-3	EndoPockets table drape	3	10
MIS1-Towel-4	EndoPockets table drape	4	10



PRO-CLIP ERCP ACCESSORY RETAINING CLIP

The autoclavable retaining clip holds your endoscopic accessories together while they are coiled up between uses.

FEATURES

- The Pro-Clip is autoclavable minimising the possibility of cross-contamination.
- Holds ERCP accessories together during the procedure.
- 10mm opening width allows accessories to be coiled up to 3 times.

Part number	Description	Pack size
MIS2-Pro-Clip	ERCP retaining clip with 10mm opening width	20

PLASTIC BILIARY STENT INSERTION KITS

Insertion kits are either available as separate items, or pre-assembled.

FEATURES

- Manufactured from ultra-smooth polyethylene material for enhanced cannulation features.
- Fluoroscopic visual marker bands to aid correct placement.
- Full range including pre-assembled individual pushers and guiding catheters for every possibility.



Part number	Diameter	Length	Type	Pack size
PR02-F5-05-220-25	8.5Fr guide catheter	220cm	Guide catheter	1
PR02-G7-08-170	8.5Fr push catheter	170cm	Push catheter	1
PR02-F5-06-220-35	10Fr guide catheter	220cm	Guide catheter	1
PR02-G7-10-170	10Fr push catheter	170cm	Push catheter	1



PLASTIC BILIARY STENT INSERTION KITS

Insertion kits are either available as separate items, or pre-assembled.

FEATURES

- Manufactured from ultra-smooth polyethylene material for enhanced cannulation features.
- Fluoroscopic visual marker bands to aid correct placement.
- Full range including pre-assembled individual pushers and guiding catheters for every possibility.

Part number	Diameter	Length	Type	Pack size
PR01-H7-05-205-32	For 8.5Fr stents	220cm	8.5Fr pre-assembled introduction kit	1
PR01-H7-06-205-35	For 10Fr stents	220cm	10Fr pre-assembled introduction kit	1



NASAL BILIARY DRAINAGE TUBE; DOUBLE PIGTAIL

The Uno-Drain ENBD tubes provide drainage of biliary ducts and/or the injection of solutions.

KIT COMPOSITION

- 1 naso-biliary tube.
- 1 nasal tube Ø 5, 7 or 10 Fr / length 50 cm.
- 1 extension containing a 3-way stopcock, length 135 cm.
- 1 removable luer-lock connector Ø 5, 7 or 10 Fr.
- 1 guidewire Ø 0,035" or 0,025" / length 400 cm.

Part number	Curve type	Ø Tube (Fr)	Length (cm)	Pack size
AS0470500	Curve 2 : double loop pigtail	5	290	1
AS0470700	Curve 2 : double loop pigtail	7	290	1
AS0471000	Curve 2 : double loop pigtail	10	290	1

GRASPING FORCEPS

The rat tooth forceps offer exceptional grasping ability of biliary stents which need removing prior to replacing.

FEATURES

- Rat tooth forceps with wide jaw opening are purpose designed for biliary stent removal.
- Can be used for straight and pigtail stents.
- Sterile packed.



Part number	Channel size	Length (cm)	Type of grasp	Pack size
FOR2-A1-23-230	2.8mm or larger	230	Rat tooth	1



CONTRAST SYRINGE

The Medical Innovations contrast syringe aids injection of thick contrast during ERCP procedures.

FEATURES

- Three ring style.
- Ergonomic design.
- 10ml capacity.

Part number	Description	Pack size
CS5801	Three ring syringe	10



CONTRAST SYRINGE

The Medical Innovations contrast syringe aids injection of thick contrast during ERCP procedures.

FEATURES

- Two finger rings.
- Ergonomic design.
- Palm plunger.
- 10ml capacity.

Part number	Description	Pack size
CC5B900	Two ring palm design	10

New product coming soon
– Affix details here

New product coming soon
– Affix details here

Products for ERCP